

EPA896Ra61 100µg
Eukaryotic Cystatin C (Cys-C)
Organism Species: *Rattus norvegicus* (Rat)
Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Eukaryotic expression

Host: 293F cell

Residues: Gly21~Ala140

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 5% Trehalose.

Original Concentration: 50µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.2

Predicted Molecular Mass: 14.9kDa

Accurate Molecular Mass: 16kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

GTSRPPRLL GAPQEADASE EGVQRALDFA
 VSEYNKGSND AYHSRAIQVV RARKQLVAGI NYLQDVEMGR TTCTKSQTNL
 TNCPFHDQPH LMRKALCSFQ IYSVPWKGTH TLTSSCKNA

[IDENTIFICATION]

ATGGAACCTCCAGGCAACCCCGCGATTGTTGGGAGCTCCGACGGAGGCAGATGCCAGGAGGGCGTGCAGCGAGCGTTGGACTTGGCCGTAAAGCGAGTACAAACAGGGCCAGCAACGATGGTACCCACAGCCGCGCCATACAGGTTGGTGGAGAGCTGGTAAGCAGCTTGTGGCTGGAAATAACTATATTTGGATGTGGAGATGGGCCGAACACTACA
 G T S R P P P R L L G A P Q E A D A S E E G V Q R A L D F A V S E Y N K G S N D A Y H S R A I Q V V R A R K Q L V A G I N Y Y L D V E M G R T T I

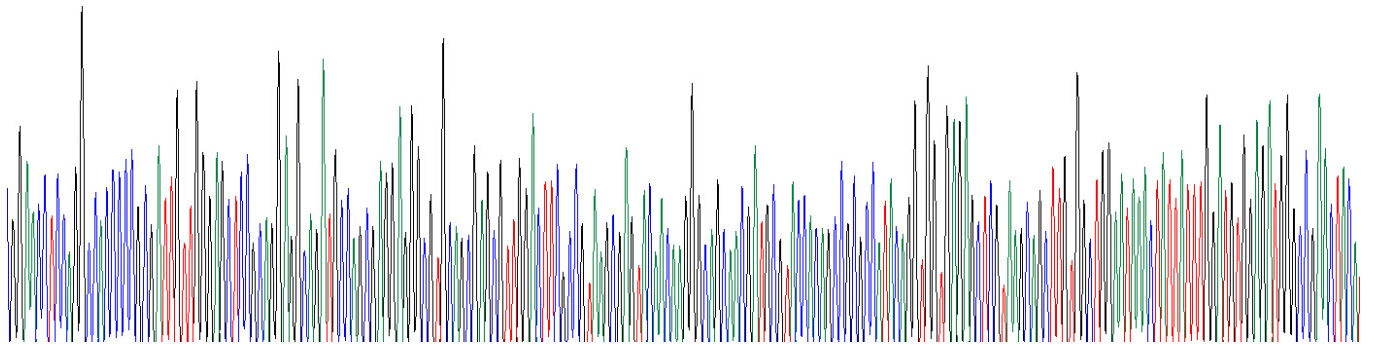


Figure . Gene Sequencing (extract)

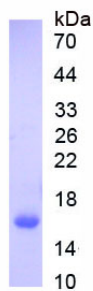


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.