

EPA033Hu61 10µg

Eukaryotic Interferon Alpha (IFN α)

Organism Species: *Homo sapiens* (Human)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Eukaryotic expression

Host: 293F cell

Residues: Cys24~Leu181

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 95%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 5% Trehalose .

Original Concentration: 350µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.9

Predicted Molecular Mass: 19.8kDa

Accurate Molecular Mass: 20kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```

CDLPETH SLDNRRTLML LAQMSRISPS
SCLMDRHDFG FPQEEFDGNQ FQKAPAI SVL HELIQQIFNL FTTKDSSAAW
DEDLLDKFCT ELYQQLNDLE ACVMQEERVG ETPLMNADSI LAVKKYFRII
TLYLTEKKYS PCAWEVVRAE IMRSLSLSTN L
    
```

[IDENTIFICATION]

```

TGTGATCTCCCTGGAGCCAGCCTGGATACAGGAGGACCTGGATCTCTCTGGACAAATGAGCGAGATCTCTCTCTCCCTCTGTGGAGGACATGCTTTGGATTTCCTGGAGGACCTTTGTTGGACAGCAGTTCGAGAGCTCCGAGCCTCTGTCTCTCATGGCTGATCGGAGATCTCAACTC
CDLPETHSLDNRRTLMLLAQMSRISPS SCLMDRHDFGFPQEEFDGNQFQKAPAI SVLHELIIQIIFNL
    
```

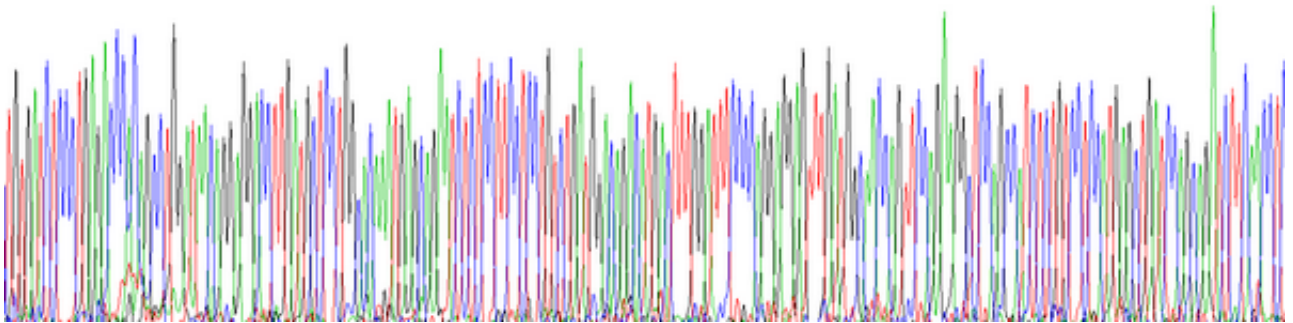


Figure . Gene Sequencing (extract)

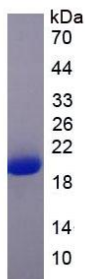


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.