

**EPA133Cp61 100µg**

**Eukaryotic Tumor Necrosis Factor Alpha (TNFa)**

**Organism Species: *Capra hircus*; *Caprine* (Goat)**

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Eukaryotic expression

**Host:** 293F cell

**Residues:** Arg79~Leu234

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted

**Purity:** > 90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 5% Trehalose.

**Original Concentration:** 300µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.5

**Predicted Molecular Mass:** 17.8kDa

**Accurate Molecular Mass:** 20kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

RS SSQASSNKPV AHVVANISAP  
 GQLRWGDSYA NALKANGVEL KDNQLVVPTD GLYLIYSQVL FRGHGCPSTP  
 LFLTHTISRI AVSYQTKVNI LSAIKSPCHR ETPEGAEAKP WYEPIYQGGV  
 FQLEKGDRLS AEINQPEYLD YAESGQVYFG IIAL

**[ IDENTIFICATION ]**

GAATT CAGGT CATCTTCTCAAGCCTCAAGCAACAAGCCGTTAGCCCAAGTGTGAGCAACATCAGGCTCCGGGCAAGCTCCGGTGGGGGACTCGTATGCCAA TGCCCTCAAGGCCAAAGGCTGGAGCTGAAAGACAACAGCTGGTGGTCCCACTGGCGGCTTTACCTCATTT

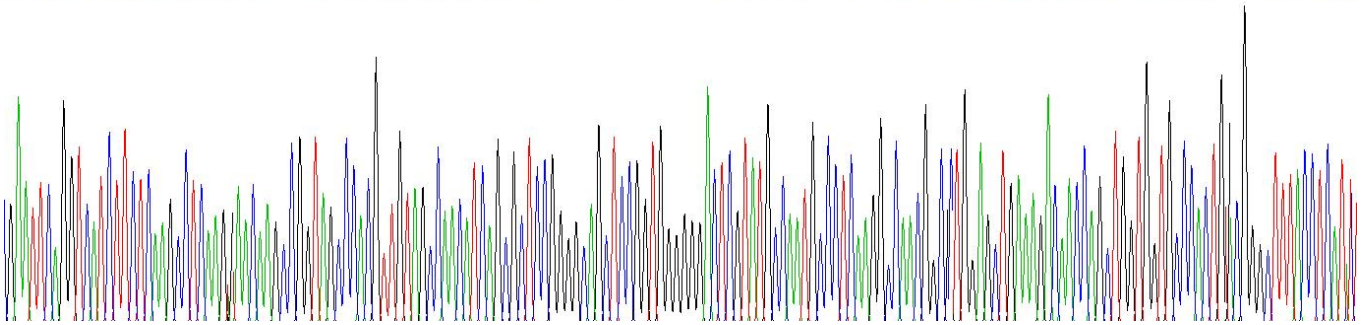


Figure . Gene Sequencing (extract)

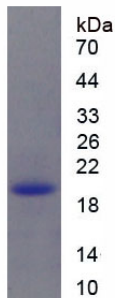


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.