

PAA183Hu01

Polyclonal Antibody to Galactosidase Alpha (GLa)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Polyclonal antibody preparation

Host: Rabbit

Purification: Antigen-specific affinity chromatography followed by Protein A affinity chromatography

Traits: Liquid

Concentration: 0.5mg/mL

UOM: 100µg(200µL)

Cross Reactivity: N/A

Applications: WB; IHC; ICC; IP; FCM

[IMMUNOGEN]

Immunogen: Recombinant GLa (Trp81~Leu429) expressed in *E.coli*

Accession No.: RPA183Hu01

[APPLICATIONS]

Western blotting: 0.01-2?g/mL

Immunohistochemistry: 5-30?g/mL

Immunocytochemistry: 5-20?g/mL

Flow cytometry:20?g/mL;

Optimal working dilutions must be determined by end user.

[FORMULATION]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 4°C for frequent use.

Aliquot and store at -20°C for 24 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined

by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[IDENTIFICATION]



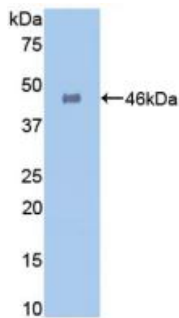



Figure. Western Blot; Sample: Recombinant GLa, Human.

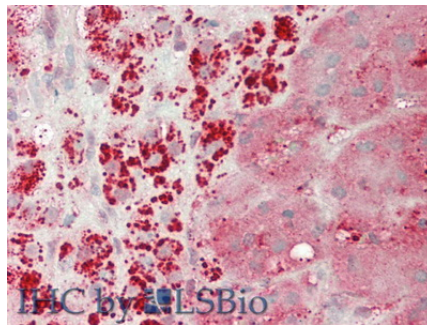

Vector Red staining on IHC-P;
Samples: Human Adrenal Gland Tissue;
Primary Ab: 15µg/ml Rabbit Anti-Human GLa Antibody
Second Ab: 2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody

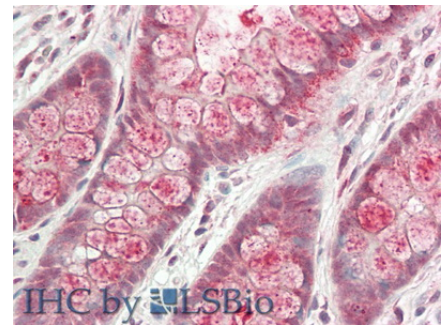

Vector Red staining on IHC-P;
Samples: Human Small intestine Tissue;
Primary Ab: 15µg/ml Rabbit Anti-Human GLa Antibody
Second Ab: 2µg/mL HRP-Linked Caprine Anti-Rabbit IgG Polyclonal Antibody

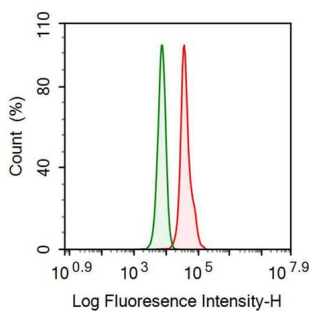

Human Hela cell was fixed with 2% paraformaldehyde (10 min) , permeabilised with 0.1% BSA-Triton X-100, then stained with 20µg/ml rabbit Anti-human GLa Polyclonal Antibody (Catalog PAA183Hu01, red histogram) or Isotype control antibody (Catalog IS067-Rb01, green histogram), followed

by 1µg/ml FITC-conjugated Anti-rabbit
IgG Secondary Antibody (Catalog
SAA544Rb18).

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.