

PAA055Po01

Polyclonal Antibody to Insulin Like Growth Factor Binding Protein 4 (IGFBP4)

Organism Species: *Sus scrofa*; Porcine (Pig)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[**PRODUCT INFORMATION**]

Immunogen: IGFBP4, Porcine

Clonality: Polyclonal

Host: Rabbit


Immunoglobulin Type: IgG

Purification: Affinity Chromatography.

Applications: WB, ICC, IHC-P, IHC-F, ELISA

Concentration: 200µg/mL

UOM: 100µg

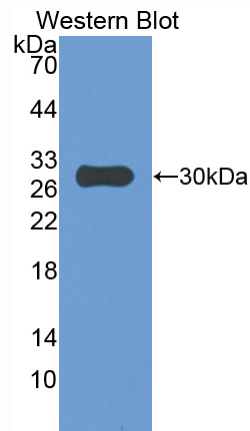

Sample: Recombinant IGFBP4, Porcine

[**IMMUNOGEN INFORMATION**]

Immunogen: Recombinant IGFBP4 (Asp23~Glu259) expressed in *E.coli*.

Accession No.: RPA055Po01

Sequence: The target protein is fused with two N-terminal Tags, His-tag and T7-tag and its sequence is listed below.

MGSSHHHHH SSGLVPRGSH MASMTGGQQM GRGSEF- DEAIHCPP CSEKLRRCR
PPVGCEELVR EPGCGCCATC ALGKGMP CGV YTPRCGSLR CYP PRGVEKP
LHTLMHGQGL CMELAEIEAI QESLQPSDKD EGDHPNNSFS PCSPQDRRCL QKHLAKIRDR
STSGGKMKVI GAPREEARPV PQGSCQSELH RALERLAASQ RRTHEDLYII PIPNCDRNGN
FHPKQCHPAL DGQRGKCWCV DRKTGVKLP GLEPKGELDC HQLADSFRE

[ANTIBODY SPECIFICITY]

The antibody is a rabbit polyclonal antibody raised against IGFBP4. It has been selected for its ability to recognize IGFBP4 in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN_3 , 50% glycerol.

[QUALITY CONTROL]

Content: The quality control contains recombinant IGFBP4 (Asp23~Glu259) disposed in loading buffer.

Usage: 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.
5uL per well when used in enhanced chemiluminescent (ECL).

Note: The quality control is specifically manufactured as the positive control. Not used for other purposes.

Loading Buffer: 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN_3 0.02%.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.