

RPJ640Hu01 50µg

Recombinant Anillin (ANLN)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Cys712~Pro1124

Tags: N-terminal His Tag

Subcellular Location: Nucleus, Cytoplasm

Purity: > 80%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose .

Original Concentration: 150µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 8.5

Predicted Molecular Mass: 50.1kDa

Accurate Molecular Mass: 50&35kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```

CQVNIKQKM QELNNEINMQ QTVIYQASQA LNCCVDEEHG
KGSLEEAEAE RLLLIATGKR TLLIDELNKL KNEGPQRKNK ASPQSEFMPS
KGSVTLSEIR LPLKADFVCS TVQKPDAANY YYLIILKAGA ENMVATPLAS
TSNSLNGDAL TFTTFTLQD VSNDFEINIE VYSLVQKKDP SGLDKKKKTS
KSKAITPKRL LTSITTKSNI HSSVMASPGG LSAVRTSNFA LVGSYTL SLS
SVGN TKFVLD KVPFLSSLEG HIYLKIKCQV NSSVEERGFL TIFEDVSGFG
AWHRRWCVLS GNCISYWYTP DDEKRKNPIG RINLANCTSR QIEPANREFC
ARRNTFELIT VRPQREDDRE TLVSQCRDTL CVTKN WLSAD TKEERDLWMQ
KLNQVLVDIR LWQPDACYKP IGKP
    
```

[IDENTIFICATION]

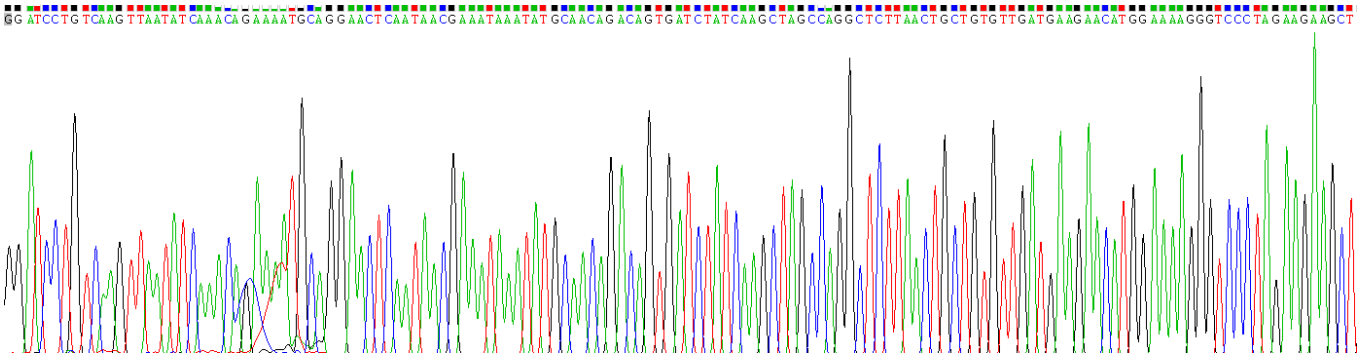


Figure . Gene Sequencing (extract)

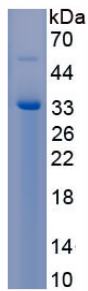


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.