

**RPJ771Mu01 10µg**

**Recombinant Cbp/p300 Interacting Transactivator, With Glu/Asp Rich Carboxy Terminal**

**Domain 1 (CITED1)**

**Organism Species: *Mus musculus* (Mouse)**

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Met1~Cys203

**Tags:** N-terminal His Tag

**Subcellular Location:** Nucleus, Cytoplasm, Chromosome

**Purity:** > 97%

**Traits:** Freeze-dried powder

**Buffer formulation:** 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% SKL, 5% Trehalose.

**Original Concentration:** 10µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 5.5

**Predicted Molecular Mass:** 24.6kDa

**Accurate Molecular Mass:** 27kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0-0.1 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.


**[ SEQUENCE ]**

MPTMSRPALD VKGGTTSKGE DANQEMNSLA YSNLGVKDRK AVTVLHYPGV  
TANGAKANGV PTSSSGSTSP IGSPTATPSS KPPSFNLHPT PHLMASMLQ  
KLNSQYQGAA ATAAAALTGA GLPGEEMPMQ NWTAPLVVG GSPGSVSPA  
GAQSPALIDS DPVDEEVLMs LVVELGLDRA NELPELWLGQ NEFDFTADFP  
SGC

**[ IDENTIFICATION ]**

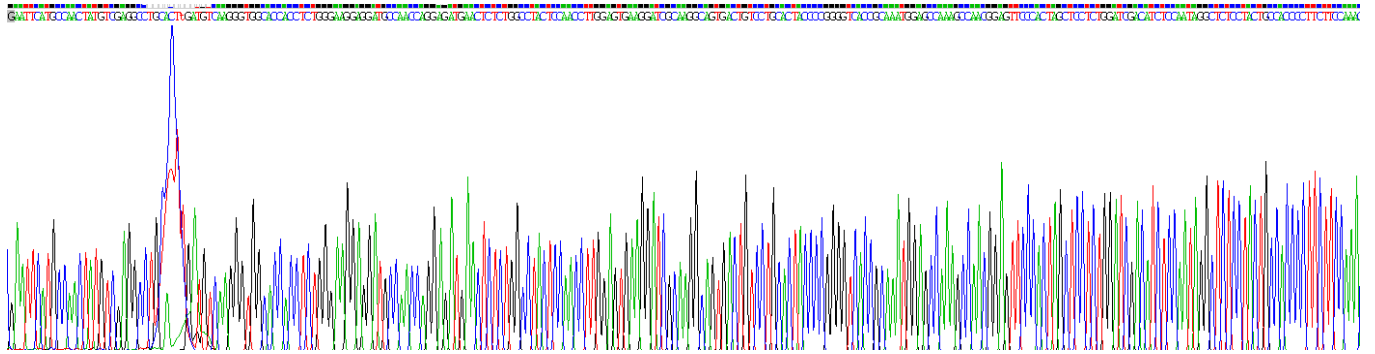

**Figure . Gene Sequencing (extract)**

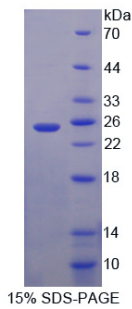

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.