

RPC474Hu03 50µg

Recombinant F-Box Protein 2 (FBXO2)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Met1~Pro296

Tags: Two N-terminal Tags, His-tag and SUMO-tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 4.2

Predicted Molecular Mass: 47.0kDa

Accurate Molecular Mass: 49kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

```
MDGDGPESV  GQPEEASPEE  QPEEASAEED  RPEDQQEEEA  AAAAAYLDEL  
PEPLLLRVLA  ALPAAELVQA  CRLVCLRWKE  LVDGAPLWLL  KCQQEGLVPE  
GGVEEERDHW  QQFYFLSKRR  RNLLRNPCGE  EDLEGWCDVE  HGGDGWRVEE  
LPGDSGVFT  HDESVKKYFA  SSFEWCRKAQ  VIDLQAEGYW  EELLDTTQPA  
IVVKDWYSGR  SDAGCLYELT  VKLLSEHENV  LAEFSSGQVA  VPQSDGGGW  
MEISHTFTDY  GPGVRFVRF  HGGQDSVYWK  GWFGARVTNS  SVWVEP
```

[IDENTIFICATION]

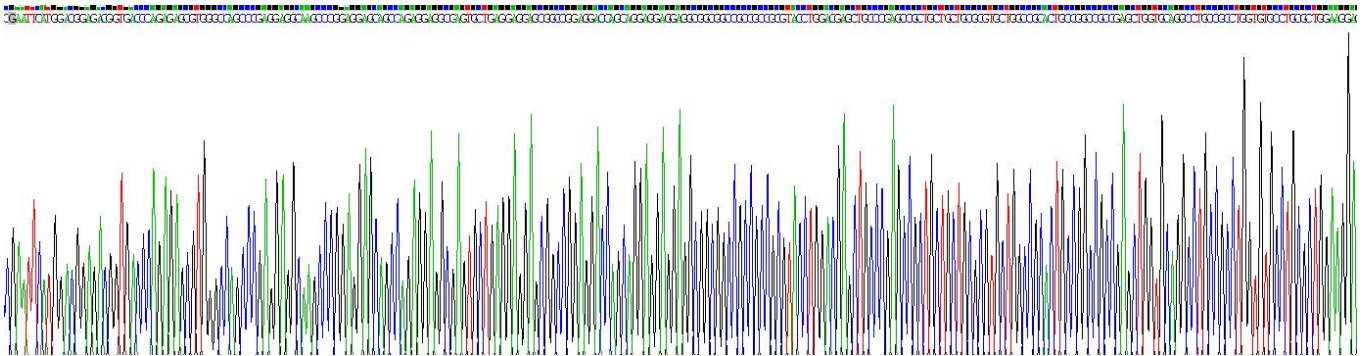


Figure . Gene Sequencing (extract)

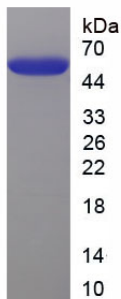


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.