

**RPA546Po01 10µg**

**Recombinant Immunoglobulin A (IgA)**

**Organism Species: *Sus scrofa*; *Porcine (Pig)***

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

12th Edition (Revised in Aug, 2016)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Ser1~Tyr341

**Tags:** N-terminal His Tag

**Subcellular Location:** Secreted

**Purity:** > 95%

**Traits:** Freeze-dried powder

**Buffer formulation:** 100mMNaHCO<sub>3</sub>, 500mMNaCl, pH8.3, containing 0.01% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 6.0

**Predicted Molecular Mass:** 40.3kDa

**Accurate Molecular Mass:** 40kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in ddH<sub>2</sub>O to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

SETSPKIFPL TLGSSEPAGY VVIACLVRDF FPSEPLVTW SPSREGVIVR  
NFPPAQAGGL YTMSSQLTLP VEQCPADQIL KCQVQHLSKS SQSVNVPCKV  
LPSDPCPQCC KPSLSLQPPA LADLLLGSNA SLTCTLSGLK KSEGVSTWQ  
PSGGKDAVQA SPTRDSCGCY SVSSILPGCA DPWVKGETFS CTAHSELKS  
ALTATITKPK VNTFRPQVHL LPPPSEELAL NELVTLTCLV RGFSPKDVLV  
RWLQGGQELP RDKYLVWESL PEPGQAIPTY AVTSVLRVDA EDWKQGDTFE  
CMVGHEALPL AFTQKTIDRL AGKPTHVNVS VVMAEAEGIC Y

**[ IDENTIFICATION ]**

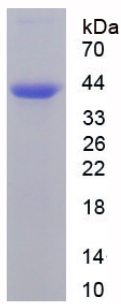


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.