

**RPF532Hu01 100 ug**

**Recombinant Trafficking Protein, Kinesin Binding 2 (TRAK2)**

**Organism Species: *Homo sapiens (Human)***

***Instruction manual***

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

---

13th Edition (Revised in Aug, 2023)

**[ PROPERTIES ]**

**Source:** Prokaryotic expression

**Host:** *E.coli*

**Residues:** Met1~Arg320

**Tags:** N-terminal His Tag

**Subcellular Location:** Mitochondrion, Cytoplasm

**Purity:** > 90%

**Traits:** Freeze-dried powder

**Buffer formulation:** PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

**Original Concentration:** 200µg/mL

**Applications:** Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

**Predicted isoelectric point:** 4.7

**Predicted Molecular Mass:** 40.7 kDa

**Accurate Molecular Mass:** 41 kDa as determined by SDS-PAGE reducing conditions.

**[ USAGE ]**

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

**[ STORAGE AND STABILITY ]**

**Storage:** Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

**[ SEQUENCE ]**

MSQSQNAIFT SPTGEENLMN SNHRDSESIT DVCSNEDLPE VELVSLLEEQ  
LPQYRLKVDI LFLYENQDWT QSPHQHQHAS DALSPVLAEE TFRYMILGTD  
RVEQMTKTYN DIDMVTHLLA ERDRDLELAA RIGQALLKRN HVLSEQNESL  
EEQLGQAFDQ VNQLQHELCK KDELLRIVSI ASEESETDSS CSTPLRFNES  
FSLSQGLLQL EMLQEKLKEL EEENMALRSK ACHIKTETVT YEEKEQQLVS  
DCVKELRETN AQMSRMTEEL SGKSDDELIRY QEELSSLLSQ IVDLQHKLKE  
HVIEKEELKL HLQASKDAQR

**[ IDENTIFICATION ]**

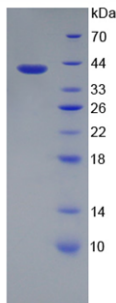


Figure. SDS-PAGE

**[ IMPORTANT NOTE ]**

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.